

Foot DETOKS device

MODEL: C01

Brand: Sinica

Operating instructions

NOT SUITABLE FOR
PREGNANT AND
Persons having received
Pacemakers

Introduction:



First you put the container with a folio (some are supplied).
The power adapter connect. The electric coil and strap attach to the device (once).
Then fill the tank supplied with water. The bracelet on the right wrist to create, so that the button is on top. Afterwards set the Elektrod rod and the feet in the water, and turn the device on (ON / OFF button).

1. Press the Start / Stop button
2. Press the MODE button. Make a fashion, you can later change.
3. The clock (Time) appears 30 minutes, you can also adjust less
4. Work voltage is 12 V
5. Current Voltage 0.9 changes etc
6. With the strap a little spoon is included.
Pour in the amount of salt in the water bath

The feet are massaged. The present in the cells toxoids (toxins) are the foot channels (more than 4000 available) released.
Only in 30 minutes the body is cleansed.

With this method, the remote toxoids, which are stored in the cell walls, purified. You should CLEANING Duch lead 14 times in succession. Each one time a day, and a total of 14 days.
If you Duch lead as a meeting 3 or 4 times a year, it is useful.

This method is used to increase physical and mental energy. The body is detoxified, the present in the cells toxins are released.

1. You achieve a safe and effective treatment with no side effects.
2. Your metabolism and immune system will be strengthened. Immune system increases the resistance of the body, increasing the risk of infections is reduced. The blood can renew itself.
3. Your memory is strengthened, bladder control is facilitated, the pH of the body is balanced. Backache, rheumatism and other pains are reduced.
4. Blood sugar, cholesterol and blood pressure can be lowered. Arthritis, gout and rheumatism can be easily detected.
5. On the body less water is accumulated, thus the weight is easily controlled. It can reduce hair loss, helps to eliminate acne and skin problems.
6. You can balance the hormonal system by negative ions, which all systems and organs are normalized, thus provides a health sign.

Quick Reference Guide

Beginning a detox therapy

1. Insert from one plastic bag (bag with some of the goods) into the bucket, so that you can easily clean the bucket. Fill the bucket with warm water up to the level of black ionizer.
2. Insert the plug of the ionizer (Elektrod) forward to the pre-forming slot (right) on the pail. Ionizer tube put into the bucket. Make sure that the ionizer coil is completely underwater.
3. Connect the power cord to the unit and to the electrical outlet.
4. Make sure you have your feet in the water and turn the unit with the red button.
5. Press the Start button and select the mode you want to run. In the first application can be eg no: Set 1. The device automatically shows 30 minutes, the duration of therapy.
6. You can use the therapy time with the - / + buttons to adjust. The maximum duration is 60 minutes. The recommended adult time is 30 minutes.
7. Check the amperage on the screen. It should not be below 2.8. You should by adding salt (eg sea salt) reach a value from 2.80 to 4.00.
8. Stir in setting the amperage the salt with the ionizer in water. Check the Amper-strength on the display.

Use little salt, so that you can use the Ioninsator for a long time.

9. If desired, you can TENS pads to the connection socket on the front panel to connect (not included)

10. If you have muscle and body aches, Fasten the hip belt to your body (hips, back, buttocks, legs). The strap to the connection socket on the front panel to connect (not included)

11. After 30 minutes, the treatment is terminated. Turn off the device with the red power button. Wash and dry the ionizer and the water tank thoroughly. The drying of the ionizer will prolong the life.

PREPARATIONS

1) Set unit on a flat surface. (The heatable hip) or the wrist strap to the device connect. In this model, only hand Schaufe.

2) Connect the power supply connector of the device with power.

3) ionizer put in the water. Attach the strap to the wrist, make sure that the place is damp on the wrist, and is on the skin.

4) To select the mode, press the MODE field, and then select the right mode.

Technical Specifications

Parameter name	Parameter Content
Supply voltage	100-240V 50/60Hz
Rated power	12,0W
current	2,80-4,00A
Time Control / Time	15-30dk.
Voltage adapter	15V/3A
Packaging	47×55×31cm

Declaration of Modfunktion

(+): Positive (-): negative

Mode 1: This mode is recommended for most people.

Mode 2: For Office Employees

Mode 3: In general, it is recommended for women.

Mode 4: Suitable for all people.

Use the (-) time (+) You can set the time. Set time is 30 minutes in the rule.

Start the cleaning with the "START" button.

5) Put your feet in the water to start the detoxification of the body.

6) With the "STOP" button to turn off the power at any time.

7) In each case, the hygiene is very important. Put the water tank with Plastitasche from (some are included).

It is recommended that you clean the electrode (ionizer), water tank and other accessories when wearing gloves. To prevent direct contact with the water after the treatment.

Colors and objects in water

Color of water	Memory place of toxoid	Basic symptoms
Weißer foam	lymphatic System	Skin allergy, irregular menstrual
White spotted	Enzyme disorder fermentation	Insomnia, constipation, flatulence
Yellow-Green	Urinary system, kidneys	Gynecological disease, prostate edema
Green	Bile, digestion system	low immune system
Brown	liver	Cigarette smoking, obesity
Orange	joints	Joint infections, rheumatism, gout
Red	Deposition in blood vessels	Narrowing of blood vessels, high blood pressure, clotting of blood
Black with spots	High accumulation of metal	Diabetes, cardiovascular disorder

Many colors and objects appear in the water during ion cleaning
When footbath a reaction between the alkalinity or acidity of the person arises while toxins and particles are released by the addition of salt in the water.

It is irrelevant whether the feet are in the water or not, these combinations cause discoloration of the water. The color of the water also depends on the geographic location. When you turn on the device without the feet, produced a reaction. The color of the water changes. This basic principle you should understand exactly.

As explained above, is the discoloration of the water by reaction between the water and all the variables that are added, depending.

This color change is dependent on the waters of specific geographical area and in the water and in the air, toxic and chemical components from.

Toxins in our bodies will look the same as the toxins that are prevalent in a particular place, but the water with their feet contains more toxic particles in comparison.

The toxoids are in the water and excreted by the body and adhesive rings rings which can usually be removed with detergent, do not stick to the water tank.

Information on the use of errors

- Do not expose to strong shocks, and leave it on the level not drop the device.
- Do not use. The unit in a hot or cold, dusty or very humid environment
- Do not expose the device to direct sunlight.
- Do not use. The product in a strong magnetic field
- Do not use to treat any chemical solutions to clean the device or rust.

Transporting the Machine

Carry your device in the original packaging. Do not drop the device to items not on, and do not let it out of hand.

For humans and the environment dangerous or damaging situation

There are no harmful situation.

Maintenance, repair or cleaning information for the consumer

The consumer does not need its own repairs carried out

Declaration of the supplied accessories



Elektro-Stap Anschluß = Electrical coil connection
Handschlaufenstecker = Wrist Strap Connector



**Elektrod-Stab/
Ionen-Stab**

Electrical coil




hand Strap



Scoops of salt

Context of the feet with the institutions:

Fußreflexzonen - rechter Fuß

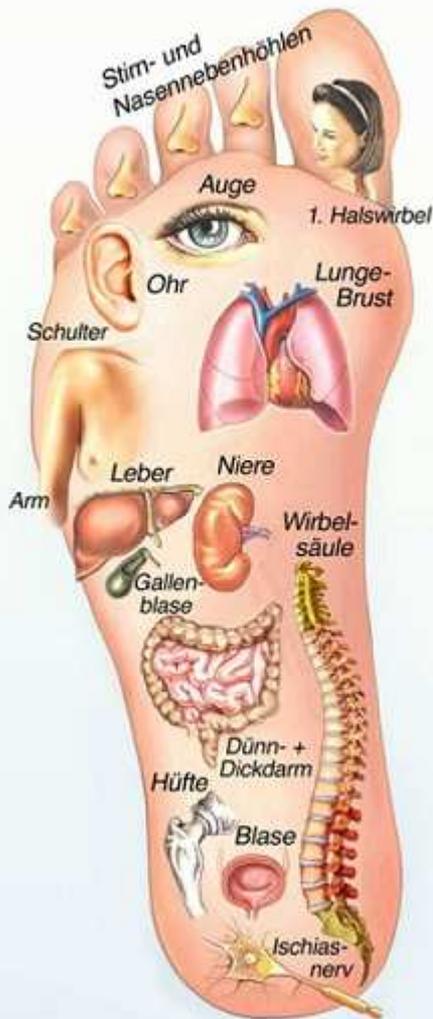

Fußreflexzonen - linker Fuß



Appearance of the water in the treatment:

